

FREEWAY *transactive* Xtra







PRISM MEDICAL UK

Moving, Handling and Bathing Solutions™


INSTALLATION SOLUTIONS

...AND A TRACKING SYSTEM FOR EVERY SITUATION

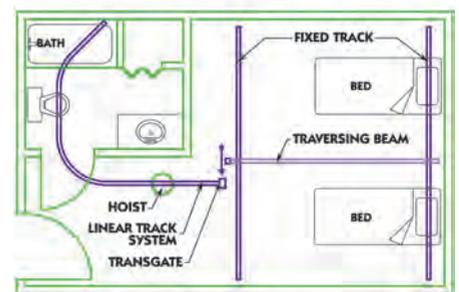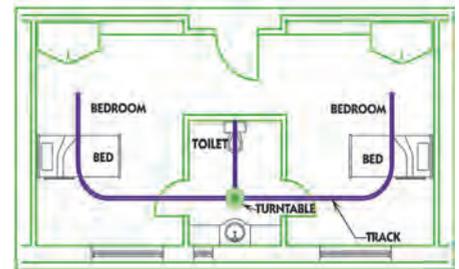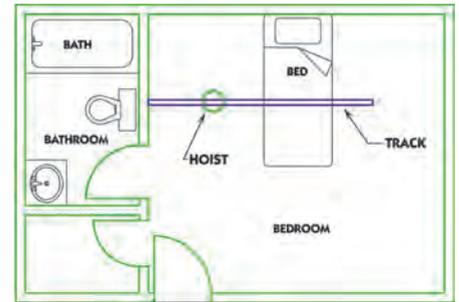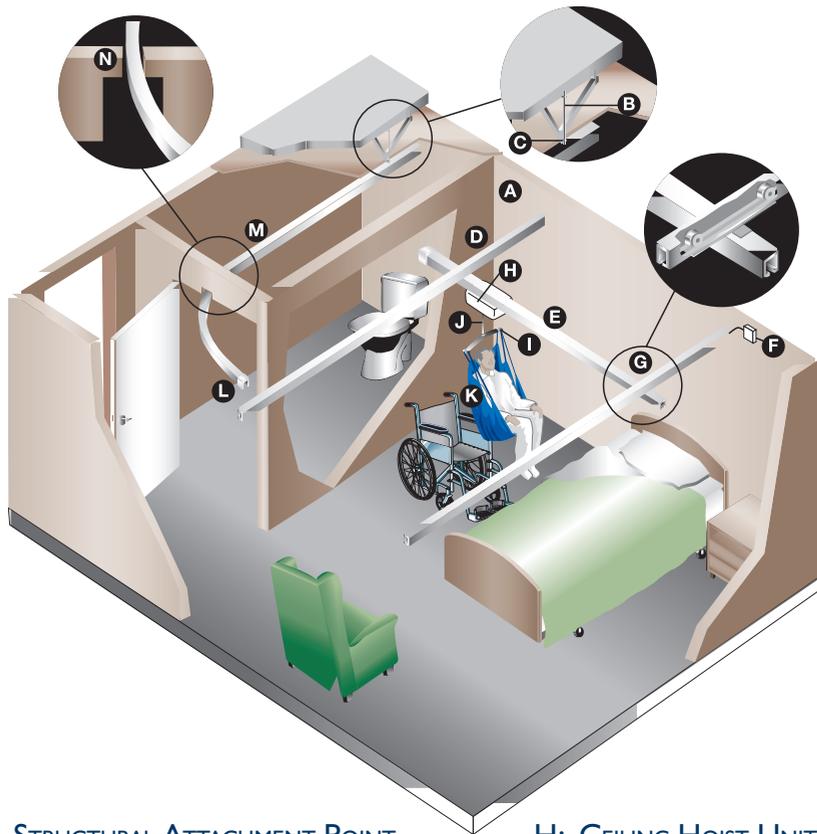

- A: STRUCTURAL ATTACHMENT POINT
- B: THREADED ROD
- C: TRACK BRACKET
- D: TRACK
- E: TRAVERSING BEAM ('H' SYSTEM)
- F: CHARGING DOCK
- G: 'H' SYSTEM

- H: CEILING HOIST UNIT
- I: CARRY BAR
- J: LIFTING TAPE
- K: SLING
- L: TRANSITION GATE
- M: CURVED TRACK
- N: MODIFIED DOOR HEADER

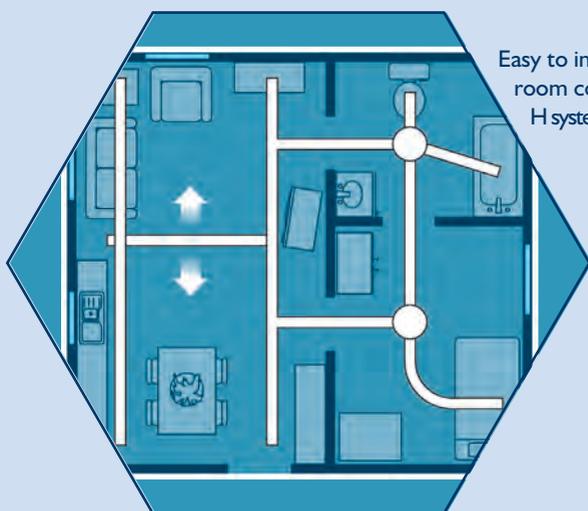

Easy to install whole room cover H systems



TRANSACTIVE XTRA

THE COMPACT HOIST WITH EXTRA FEATURES...

The all new Transactive Xtra, a product of world class design and state of the art features providing safe, reliable and effective hoisting for carers and users alike. Together with the extensive range of track configurations, the Transactive Xtra can provide the perfect solution to your hoisting requirements.

- **Slimline and with great lifting capacity**
- **Low profile hoist unit is unobtrusive and provides greater lifting height**
- **All lifting and traversing operations have soft start/stop to ensure smooth transfers for the user**
- **A host of safety features**
- **Smooth running tape with large lifting range**
- **State of the art display with lift counter**

Slimline low profile padded carry bar has six sling loop attachment points as standard (four hook as option) different width carry bars available

Ergonomically designed handset with hold to run clearly marked buttons and directional arrows requires minimum pressure to operate. Pneumatically controlled, it is therefore safe for use in wet conditions.

Fire retardant cover completely encapsulates hoist

Neat in track transformer IP 67 rated

Coloured Direction indicators on hoist unit match hand control to ensure correct operation first time


Emergency stop and electric emergency window

**130kgs, 160kgs, 200kgs, 270kgs
MODELS AVAILABLE**



Display provides essential information for the user carer and engineer.

- Lift counter
- Displays operation being carried out (UP/ DOWN etc)
- Charge level of battery • Low battery indicator (visual and audial)
- On charge indicator • Return to charge (Option)
- Diagnostic feature for transactive engineer to aid maintenance and repair

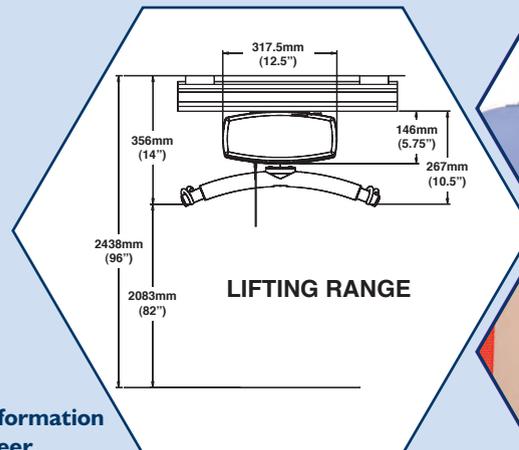


Manual window and raising.

This feature is unique to Freeway hoists and allows the hoist to continue to be used until the engineer arrives.

This provides peace of mind and security to carer and user.

TRACK AND ACCESSORIES

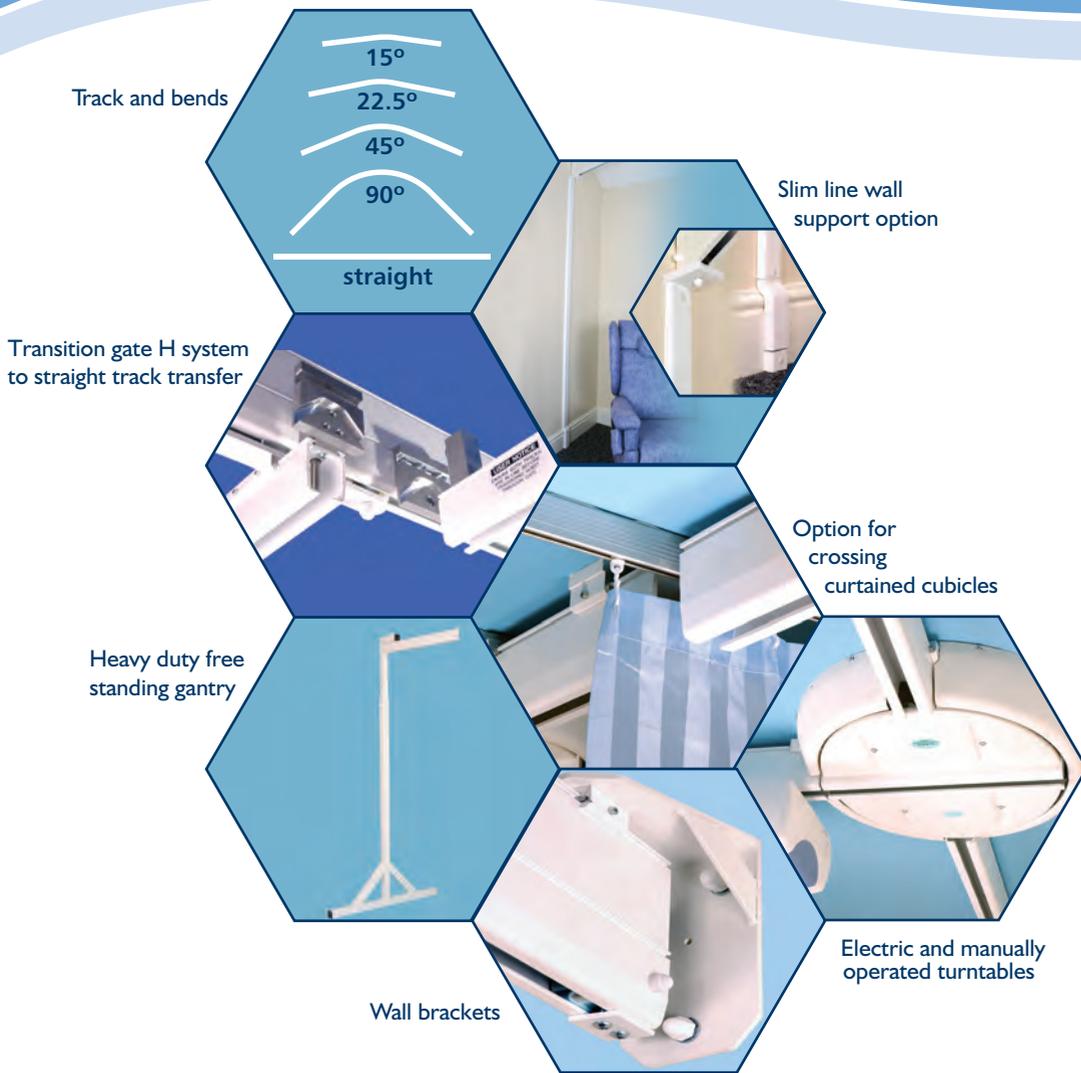

PRISM MEDICAL UK

Prism Medical UK Divisions:

Community Care

Prism Medical UK – Westholme Division (Northern Region)

Unit 4 Jubilee Business Park,
Jubilee Way, Grange Moor, Wakefield,
West Yorkshire WF4 4TD UK
Tel: +44 (0) 844 980 2260
Fax: +44 (0) 844 980 2261
E-mail: info@westholme.co.uk
www.prismmedical.co.uk

Prism Medical UK – Westholme Division (Southern Region)

Unit 7, Saracen Industrial Estate,
Mark Road, Hemel Hempstead,
Hertfordshire HP2 7BJ UK
Tel: +44 (0) 844 980 2263
Fax: +44 (0) 844 980 2264
E-mail: info@westholme.co.uk
www.prismmedical.co.uk

Prism Medical UK – Liftech Division & Test Valley Division

Unit 42 South Hampshire Industrial Park, Totton,
Southampton, Hampshire SO40 3SA UK
www.prismmedical.co.uk

Liftech Division

Tel: +44 (0) 844 980 2272
Fax: +44 (0) 844 980 2273
E-mail: sales@liftechsystems.co.uk

Test Valley Division

Tel: +44 (0) 844 980 2269
Fax: +44 (0) 844 980 2270
E-mail: admin@testvalleymobility.co.uk

Prism Medical UK – Carefree Division

13 - 15 Pennington Street, Hindley,
Wigan, Lancashire WN2 3AH UK
Tel: +44 (0) 844 980 2266
Fax: +44 (0) 844 980 2267
E-mail: info@carefreelifting.co.uk
www.prismmedical.co.uk

Prism Medical UK - Movement2 Division

Unit 8 The Warren,
East Goscote Industrial Estate,
East Goscote, Leicestershire LE7 3XA
Tel: +44 (0) 844 980 2281
Fax: +44 (0) 844 980 2282
E-mail: mail@movement2
www.movement2.co.uk

Long Term Care

Prism Medical UK – Long Term Care

Unit 4 Jubilee Business Park,
Jubilee Way, Grange Moor, Wakefield,
West Yorkshire WF4 4TD UK
Tel: +44 (0) 844 980 2260
Fax: +44 (0) 844 980 2261
E-mail: carcinfo@westholme.co.uk
www.prismmedical.co.uk

Prism Medical UK – Saluss Division

Unit 1, Tir Llwyd Industrial Estate, St Asaph Avenue,
Kinmel Bay, Rhyl, Conwy LL18 5JA UK
Tel: +44 (0) 844 980 2290
Fax: +44 (0) 844 980 2291
E-mail: helpdesk@saluss.co.uk
www.prismmedical.co.uk

Acute Hospital Care

Prism Medical UK – Healthcare Division

Unit 4, Crispin Court, Newbiggin Lane, Newcastle Upon
Tyne, Tyne & Wear NE5 1BF UK
Tel: +44 (0) 844 980 2278
Fax: +44 (0) 844 980 2279
E-mail: healthcareinfo@prismmedical.co.uk
www.prismmedical.co.uk

Prism Medical UK – APLS Division

Unit 4, Crispin Court, Newbiggin Lane,
Newcastle Upon Tyne, Tyne & Wear NE5 1BF UK
Tel: +44 (0) 844 980 2275
Fax: +44 (0) 844 980 2276
E-mail: info@assured-pls.com
www.prismmedical.co.uk

SPECIFICATION

MAX SAFE WORKING LOAD: 130KGS,
160KGS, 200KGS, 270KGS.

APPROX NUMBER OF LIFTS BEFORE
RECHARGE: WITH 85 KG LOAD APPROX 65

CHARGING TIME: 8HOURS MAX

UNIT WEIGHT INC BATTERIES:
MANUAL: 9.7KGS
DRIVE: 10.4KGS

LIFTING SPEED:
LOADED 0.15M/S; UNLOADED 0.25M/S

BATTERIES: 2 X 12V 5AMP
BATTERY PROTECTION: SLEEP MODE
VOLTAGE 24VDC

TRANSFORMER : IP 67 RATED
30V AC OUTPUT
230V AC INPUT
50/60MHZ



FULL 2 YEAR WARRANTY (excludes batteries)

The Freeway TransactiveXtra is manufactured in the UK to meet and exceed all relevant design, safety and quality standards.

✓ BS EN ISO 10535:2006 Compliant 'Hoists for the transfer of disabled persons – Requirements and test methods'



PRISM MEDICAL UK



Unit 1, Tir Llwyd Industrial Estate, St Asaph Avenue, Kinmel Bay, Rhyl, Conwy, LL18 5JA UK

Tel: +44 (0) 845 409 2000 E-mail: info@prismmedical.co.uk www.prismmedical.co.uk